



How to use **Awesome Floss**[®] dental floss*


Front Teeth

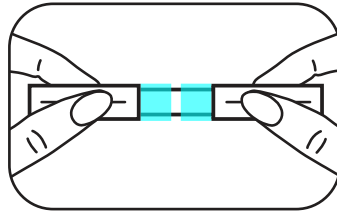


1 Grab handles only.

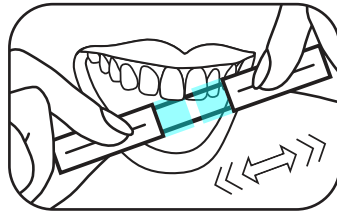


2 Glide one strand back and forth between teeth.
(Top strand is for upper teeth and bottom strand is for lower teeth.)

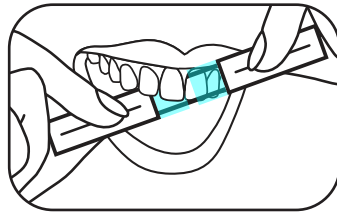


3 Ease that strand gently under gums.

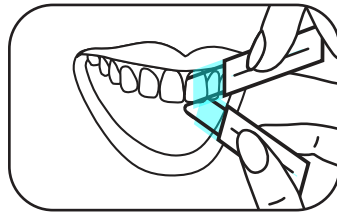

4 Hug that strand tightly against one tooth.

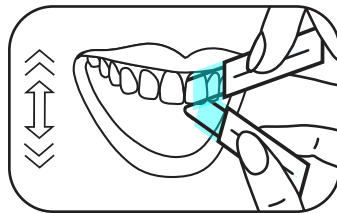


5 Gently move hugging strand up and down several times to clean properly.

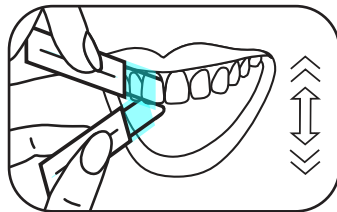

6 Repeat Steps 4 and 5 on adjacent tooth.

7 Repeat Steps 2 thru 6 for every pair of adjacent teeth.

Back Teeth

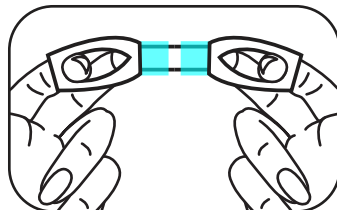

8 For back teeth: Insert forefingers in slits of handles. Press tips of middle fingers against handles.